


Pressure Vessel


Closure Style Comparison

Parker Autoclave Engineers' offers three basic vessel closure styles. Select the style that meets your requirements for ease of opening and closing, seal compatibility, as well as pressure and temperature performance.

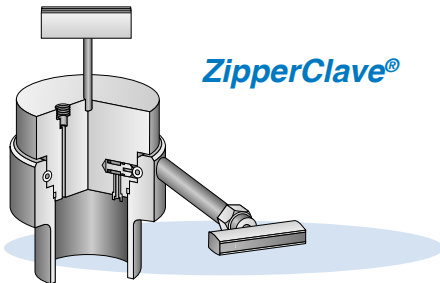


ZipperClave®

The ZipperClave® offers the quickest opening autoclave available. No bolts to torque, nor clamps or rings. Inserting a single spring section closes the vessel. The ZipperClave® uses an elastomeric O-ring for its pressure seal.

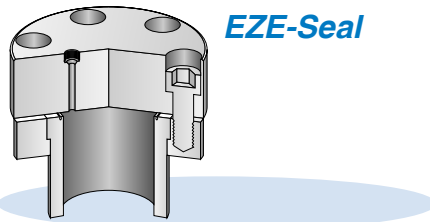

EZE-Seal

The EZE-Seal provides the ability to operate at high temperature and moderate pressure. The “loose flange” allows for easy interchange of vessel bodies. The seal design requires low bolting torque. The EZE-Seal may be equipped with either metal or elastomeric pressure seal.

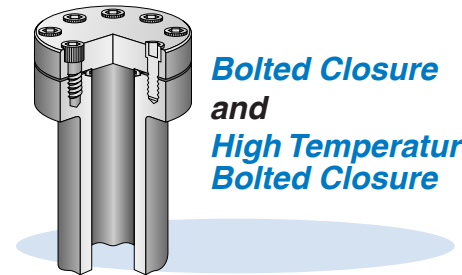

Bolted Closure and High Temperature Bolted Closure

The Bolted Closure [BC] and High Temperature Bolted Closure [HTBC] offer a high pressure capacity at moderate and high temperatures. The Bolted Closure uses multiple gasket materials as the High Temperature Bolted Closure utilizes the Gasche metal gasket for its pressure seal.

MAWP: 2,200 psi @ 450° F
(151 bar @ 232° C)

Sizes: 500 ml
1,000 ml
2,000 ml
4,000 ml

Seal Materials:
Buna N
Ethylene-Propylene
PTFE
Viton®
Silicone
Kalrez®

Vessel Materials:
ANSI 316 Stainless Steel
Hastelloy® C-276

Notes: MAWP based upon Viton® seals, temperature and pressure ratings will vary depending upon seal material selected. See Bulletin “PV-ZIP” for details of vessel performance. ASME Code Stamp not available for the Zipper Closure. Alternative materials are available, please consult factory.

MAWP: 3,300 psi @ 850° F
(227 bar @ 454° C)

Sizes: 100 ml
300 ml
500 ml
1,000 ml
2,000 ml
4,000 ml

Seal Materials:
Metal Double Delta
Buna N
Ethylene-Propylene
PTFE
Viton®
Silicone
Kalrez®

Vessel Materials:
ANSI 316 Stainless Steel
Hastelloy® C-276

Notes: MAWP based upon metal seals, temperature and pressure ratings will vary depending upon seal material selected. See Bulletin “PV-EZE” for details of vessel performance. Alternative materials are available, please consult factory.

MAWP: 5,500 psi @ 650° F [BC]
(379 bar @ 343° C)
5,000 psi @ 950° F [HTBC]
(345 bar @ 510° C)

Sizes: 100 ml, 300 ml, 500 ml
1,000 ml
2,000 ml
4,000 ml, 8,000 ml

Seal Materials: [BC - All Materials]
Metal Gasket [HTBC only]
Buna N
Ethylene-Propylene
PTFE
Viton®
Silicone
Kalrez®

Vessel Materials:
ANSI 316 Stainless Steel
Hastelloy® C-276

Notes: MAWP based upon metal seals, temperature and pressure ratings will vary depending upon seal material selected. See Bulletin “PV-BC” for details of vessel performance. Alternative materials are available, please consult factory.

WARNING

FAILURE, IMPROPER SELECTION OR IMPROPER USE OF THE PRODUCTS AND/OR SYSTEMS DESCRIBED HEREIN OR RELATED ITEMS CAN CAUSE DEATH, PERSONAL INJURY AND PROPERTY DAMAGE.

This document and other information from Parker Hannifin Corporation, its subsidiaries and authorized distributors provide product and/or system options for further investigation by users having technical expertise. It is important that you analyze all aspects of your application and review the information concerning the product or system in the current product catalog. Due to the variety of operating conditions and applications for these products or systems, the user, through its own analysis and testing, is solely responsible for making the final selection of the products and systems and assuring that all performance, safety and warning requirements of the application are met. The products described herein, including without limitation, product features, specifications, designs, availability and pricing, are subject to change by Parker Hannifin Corporation and its subsidiaries at any time without notice.

Offer of Sale

The items described in this document are available for sale by Parker Hannifin Corporation, its subsidiaries or its authorized distributors. Any sale contract entered by Parker will be governed by the provisions stated in Parker's standard terms and conditions of sale (copy available upon request).



Instrumentation Products Division
Autoclave Engineers Operation
8325 Hessinger Drive
Erie, Pennsylvania 16509-4679 USA
PH: 814-860-5700 FAX: 814-860-5811
www.autoclaveengineers.com

Parker Hannifin Manufacturing Ltd.
Instrumentation Products Division, Europe
Industrial Estate Whitemill
Wexford, Republic of Ireland
PH: 353 53 914 1566
FAX: 353 53 914 1582

